

# LOR Balloon Indicator Syringe

## Features:

- **Visibility** : LOR Indicator Syringe makes the puncture process visualized.
- **Safety** : The pressure in balloon is low and stable.
- **Accuracy and reliability** : The combination of visible indication and experience improves the accuracy and reliability.
- **Simple operation** : Reduce the actions of pushing the syringe repeatedly.
- **Low resistance** : The use of rubber plug and its dynamic fit with the barrel make the injection smooth.
- **Clinical teaching** : The visibility can make the teacher and student judge the depth of the needle tip



“A visible “brake signal” , Balloon, let us know whether the epidural needle arrives at the epidural space”

**TUORen**

• Tuoren Medical Device India PVT. LTD.

#1712, 17th Floor Tower-B, Emaar Digital Greens, Sector - 61, Golf Course Ext. Road , Gurgaon, 122011.

• +91 124 4003250

• info@tuoren.in

• [www.tuoren.in](http://www.tuoren.in)

**TUORen**

*Redefining Healthcare Innovation*

**EpiLOR™**

**Anesthesia Kit**



With Special LOR Balloon Indicator Syringe



• [www.tuoren.in](http://www.tuoren.in)



- 1- Spinal Needle
- 2- Epidural Needle
- 3- Hypodermic Needle 18G
- 4- Hypodermic Needle 22G
- 5- Hypodermic Needle 25G
- 6- Hypodermic Syringe
- 7- LOR Balloon Indicator Syringe
- 8- Catheter Fixing Device
- 9- Snap Connector
- 10- Anaesthesia Catheter
- 11- Drug Filter


**EpiLOR™ Pro Specifications:**

Type	Specifications	OD (G)	Length (mm)	Product Code	
EpiLOR™ Pro Indicator	Combined Kit with Indicator	16G/25G	80 mm	80mm*113mm	ICML1625
EpiLOR Pro™ Indicator	Combined Kit with Indicator	18G/27G	80 mm	80mm*113mm	ICML1827
EpiLOR Pro™	Combined Kit without Indicator	16G/25G	80 mm	80mm*113mm	ICMN1625
EpiLOR Pro™	Combined Kit without Indicator	18G/27G	80 mm	80mm*113mm	ICMN1827

**Features:**

- **LOR Balloon indicator syringe** makes the puncture process visualized, improves success rate and safety of puncture.
- **Anti-injury anesthesia catheter** has strong tensile property, blue soft tip reduces injury during placement.
- **The epidural needle** deals with special process, has clear puncture feeling and smooth insertion of anesthesia catheter.
- **Spinal needle with Locking type** can be fixed on the epidural needle by a twisting-lock system to prevent the spinal needle from moving during injection. Digital indication on the epidural needle make it convenient for judging the puncture depth
- **0.2µm liquid filter**, precision filtration.
- **Snap Connector**, Ergonomic design single click closed & connection with catheter without leakage.
- **Catheter Fixing Device**, is used to help stop unintended catheter migration, Reduces kinking of the catheter leaving the body thanks to solid puncture site cover.


- **Atraumatic design catheter** , Kink – resistant, smooth tip formation to reduce insertion-related trauma.
- **Snap Connector**, Ergonomic design single click closed & connection with catheter without leakage.
- **Catheter Fixing Device**, patient discomfort is reduced because of the large foam and medical grade adhesive pad.
- **Tuohy Needle** is designed to offer precise and trustworthy access to the epidural space.

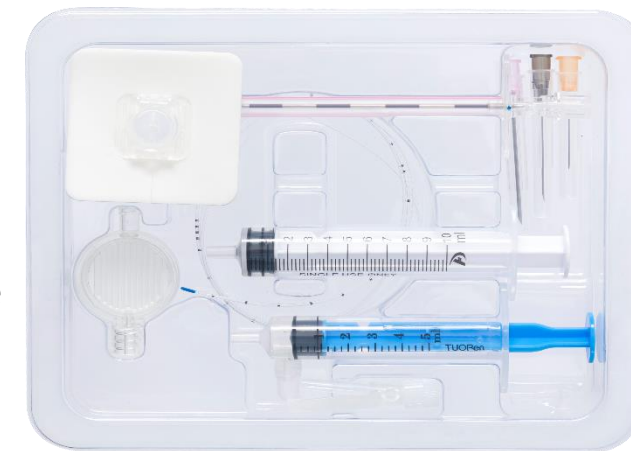

**EpiLOR™ Specifications:**

Type	Specifications	OD (G)	Length (mm)	Product Code
EpiLOR™ Indicator	Combined Kit with Indicator	16G	80 mm	IEML16080
EpiLOR™ Indicator	Combined Kit with Indicator	18G	80 mm	IEML18080
EpiLOR™	Combined Kit without Indicator	16G	80 mm	IEMN16080
EpiLOR™	Combined Kit without Indicator	18G	80 mm	IEMN18080

**Spinal Kit**

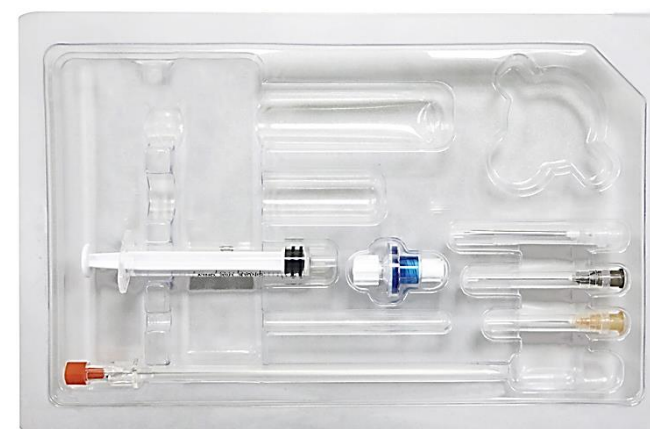

- **Pencil-point spinal needle** dealt with special process, makes puncture point heal quickly and automatically and reduces postoperative headache rate.
- It reduces the occurrence of **Postdural Puncture Headache (PDPH)** .
- **0.2µm liquid filter**, precision filtration.

**Specifications:**

OD (G)	Length (mm)	Product Code
20G	90 mm	STM2090
22G	90 mm	STM2290
25G	90 mm	STM2590
	113mm	STM25113